

Maxi Personal amplifier

© Bellman & Symfon, 2013

bellman.com

Maxi

Personal amplifier

Digital | User friendly | Robust



Unbeatable for **group 2** users

Group 1

Hard of hearing
No other
impairment



Group 2

Hard of hearing

+

Low vision

Limited dexterity

Tremors

Listening **situations**



**One-to-one
conversations**



**Theater, church
concert hall etc.**



**Listening to the
TV or radio**



**Shop / pharmacy
or doctor**



**Family dinner
restaurant / café**




In the car

Hear **even** in busy surroundings

Listen at the volume of your choice without others having to raise their voices.

Maxi will amplify speech and remove disturbing background noise.

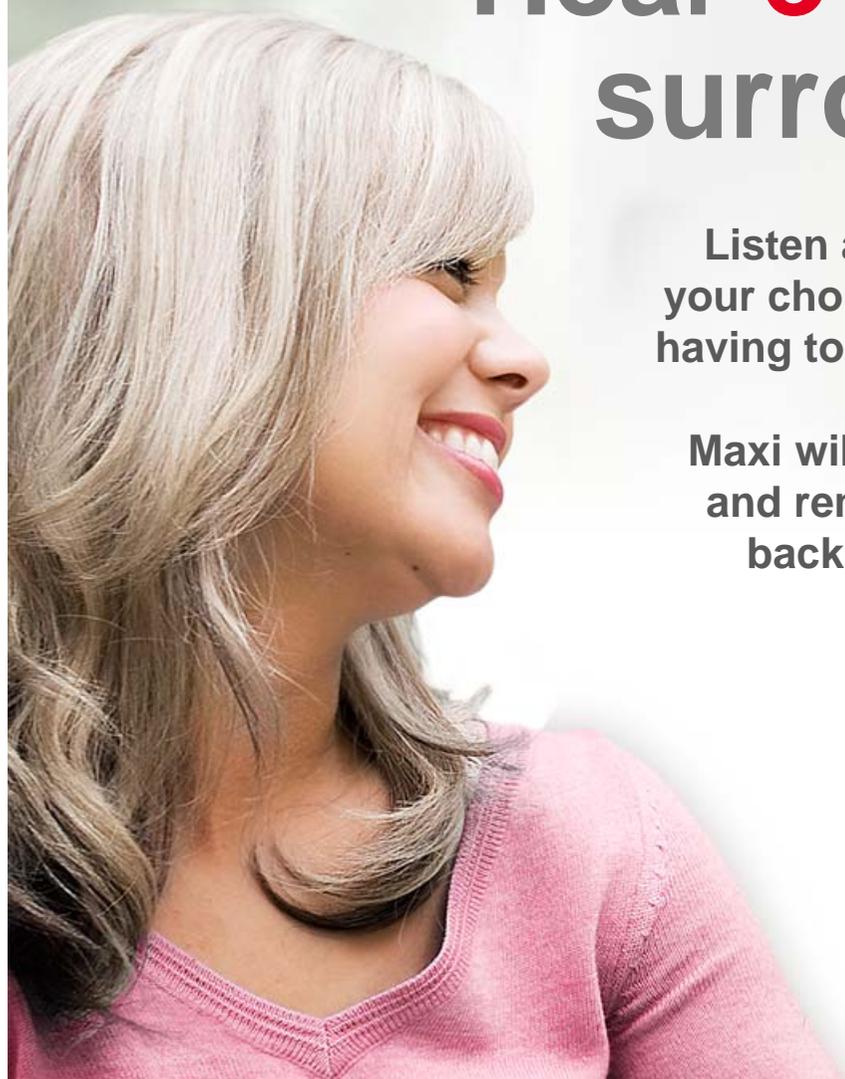

Enjoy radio and TV in full **stereo**



Connect Maxi to the TV or radio and enjoy great sound in stereo while sparing your family's ears.

Press the **push-to-listen** button when you want to talk to the person next to you.



Maxi | In theaters



Get **more** out of the show

Many theaters, concert halls, and churches are fitted with hearing loop systems that Mino picks up with its built-in **Telecoil**.

Maxi | Shop, pharmacy or doctor



Never any confusion

When the advisor's voice blends in with the background it is easy to miss important information.

Place Maxi on the counter to hear the advisor loud and clear, without her having to raise the voice.



Hear above the **traffic** noise

A car is one of the hardest sound environments to conduct a conversation in. Besides the traffic noise, the passengers are normally facing away from you.


Attach a **lapel microphone** to your travel companion's sleeve.



Easy to use

Up and running in minutes

- No complicated menus to learn
- Just connect the headphones or neck loop and adjust the volume to your liking
- One button - one function

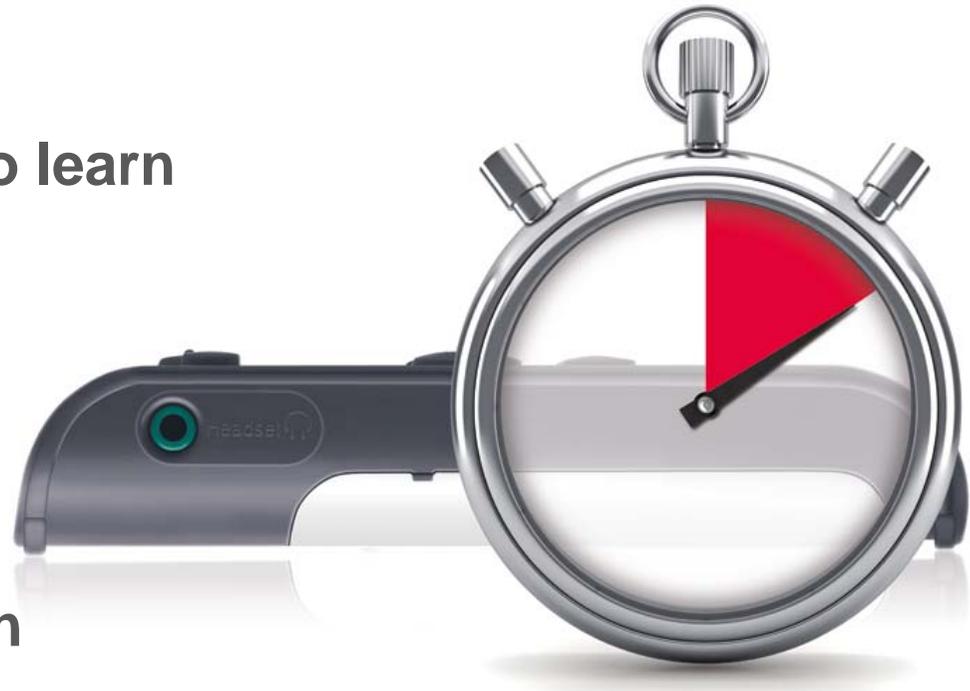

Great ergonomics

High quality and robust

- **Very robust form factor**
- **Large tactile buttons**
Unique soft grip casing
High contrast indicators
Lockable buttons
- **A great choice for those with low vision or limited dexterity**



Clear and distinct sound

Listen without effort

- Digital signal processing
- 10 band noise cancellation
- 10 band dynamic compression
- Adaptive feedback cancellation
- Sound with natural colour



Using **Maxi** with hearing aids



1. On the table

Mino captures the sound, filters out background noise and clarifies speech.

2. Around the neck

The neck loop transmits the speech wirelessly to the hearing aid or cochlear implant.


3. In the ear

The hearing aid or cochlear implant refines the sound further before it reaches the ear.

Connections



Buttons & controls

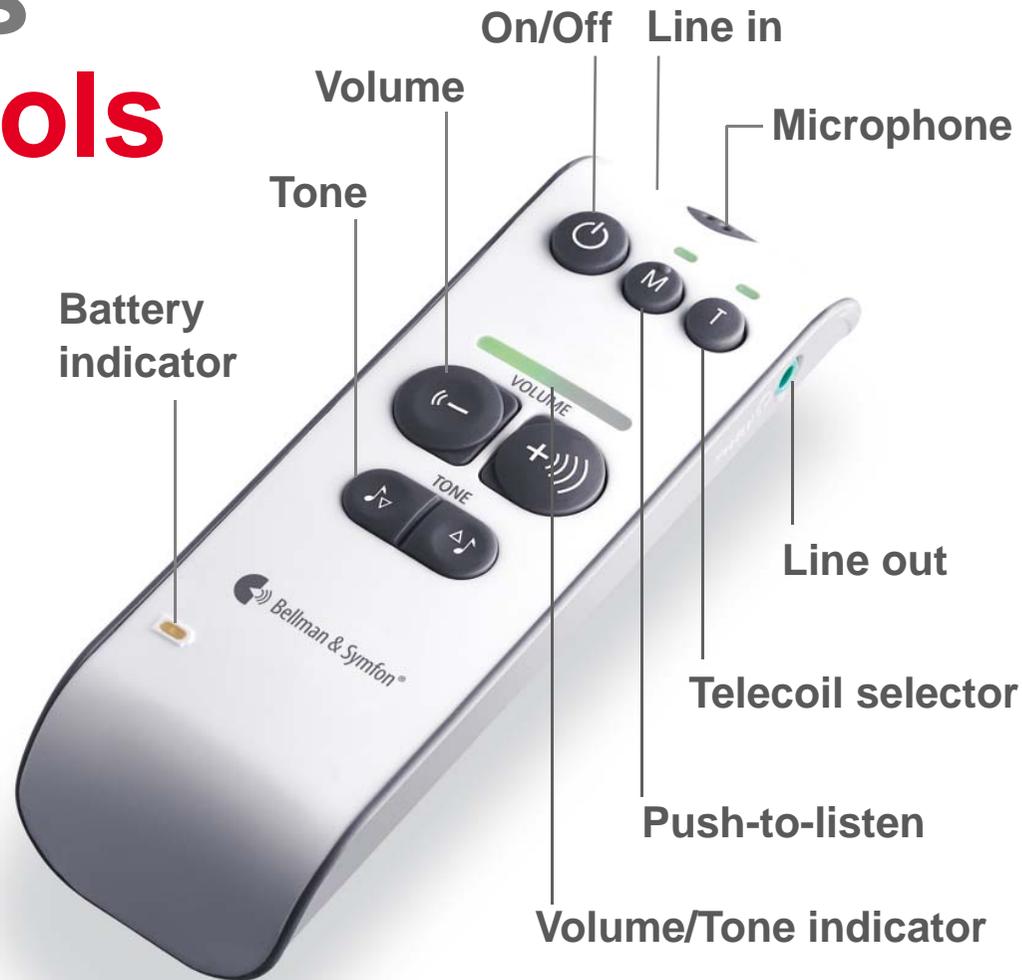

Advanced **settings** menu



Balance control

Compensates for example for hearing loss in one ear.



Basic attenuation

Trim the attenuation to a specific hearing accessory.



Disable Telecoil

Makes it easier for the end users to handle Maxi.



Tone control blocking

Makes it easier for the end user to handle Maxi.



Tone control preset

Select preset tone. Used with Tone control blocking.

Accessories



Tie clip microphone



Cable kit 5m, 16"



Neck loop



Travel case



Head phones



Ear phones



Stetoclips

Powerful **Neck** loop

- High output power
- Fabric tube for good comfort
- Safety break-away neck strap
- New pliable cable material
- Tie-clip for optimum positioning
- Soft spray-coated housing



Addressing real life needs (recap)



Digital signal processing

→ Noise free, natural sounding speech



Robust and ergonomic

→ Large buttons, soft grip, high contrast indicators



Intuitive to use

→ No complicated menus to learn

Addressing real life needs (cont.)



Built-in Telecoil

→ For venues fitted with a loop system



Flexible

→ Connects to the TV, radio, stereo system etc.



Quality of Sweden

→ Comes with a two year worldwide warranty

